




PRACTICAL INNOVATION

Ziath

DataPaq™ Multi-Rack Scanner

Fast, simple to install and with a 2-year guarantee

The Ziath DataPaq™ Multi-Rack Scanner scans SBS and cryoracks. With the choice of five different masks, nearly any combination of racks is possible and up to four racks can be scanned at the same time.

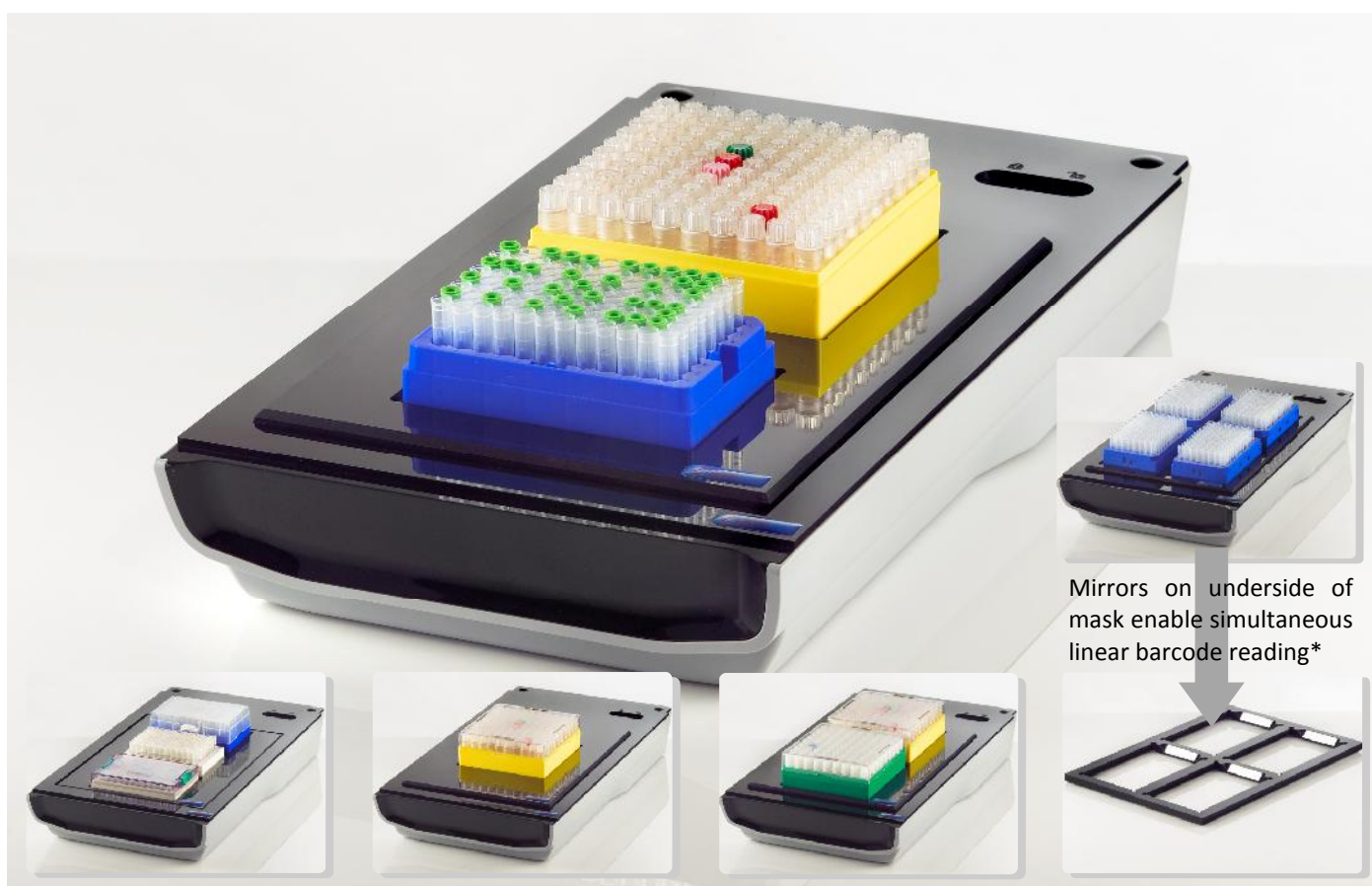

Features

- Compatible with SBS and cryoracks
- Scan SBS & cryoracks at the same time
- Comes as standard with one mask configuration
- Compatible with Matrix, ABgene, REMP, Micronic, FluidX/FCR tubes, Wheaton, Greiner and others
- Scan-time is from 5 seconds per rack
- Export data to Excel, text, XML, images
- Operating systems: Windows 7, XP, Vista
- Connects with most external databases
- Dimensions: 271 x 475 x 73mm (WxDxH)
- Durable construction - 2 year guarantee
- Free software upgrades for life
- Ziath's expertise and support come as standard

**Linear barcode must be placed on lower edge of rack skirt*




PRACTICAL INNOVATION

DataPaq™ Multi-Rack Scanner

Unique Configuration Options

The DataPaq™ Multi-Rack Scanner offers the unique ability to adapt to the configuration required. It comes with one mask of choice as standard and other masks can be added as required. See the ordering information, below, for a full list of mask configurations available.

The 4 x SBS configuration automatically reads the linear barcode on the rack skirt, by the clever use of mirrors on the underside of the configuration mask.

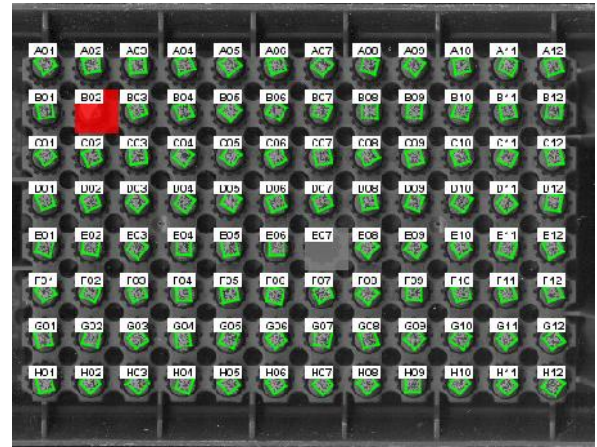

DataPaq™ software gives intelligent scan information.

Easy to Set-Up

After a one-time automatic calibration, the scanner is ready to read Matrix, ABgene, FluidX/FCR, REMP, Wheaton & Greiner tubes. The scan time is from five seconds per rack, even if a combination of different manufacturers' racks are scanned at the same time.

Our DataPaq™ software makes it simple to export to Excel, XML, JSON or text and scanned images can also be saved.

Intelligent Scan Information

The software displays the results of the scan and can identify the difference between a correctly scanned tube (green), a damaged tube (red) and missing tube (grey).

Cryoprotection™ Option

The occurrence of condensation whilst scanning tubes straight from cryogenic storage is well known. Ziath's cryoprotective coating eliminates the occurrence and ensures uninterrupted scanning. As Ziath's method for eliminating condensation uses no extra electrical components, reliability is increased and costs are significantly decreased.

Ordering Information

Code	Description	Quantity
ZTS-A4	DataPaq™ Multi-Rack Scanner	1
ZTS-A4-80	DataPaq™ Multi-Rack Scanner with Cryoprotection™	1
ZTS-MSB3	Mask for 3 x SBS Racks	1
ZTS-MCRYO2	Mask for 2 x Cryoracks	1
ZTS-MSBCR	Mask for 1 x SBS Rack and 1 x Cryorack	1
ZTS-MCRYO1	Mask for 1 x Cryorack	1
ZTS-MSB4	Mask for 4 x SBS Racks	1